
**Dr. Mayer**  
LIFE & HEALTH

## Rechargeable Sonic Power Toothbrush

Operating Instructions

Model: GTS2060

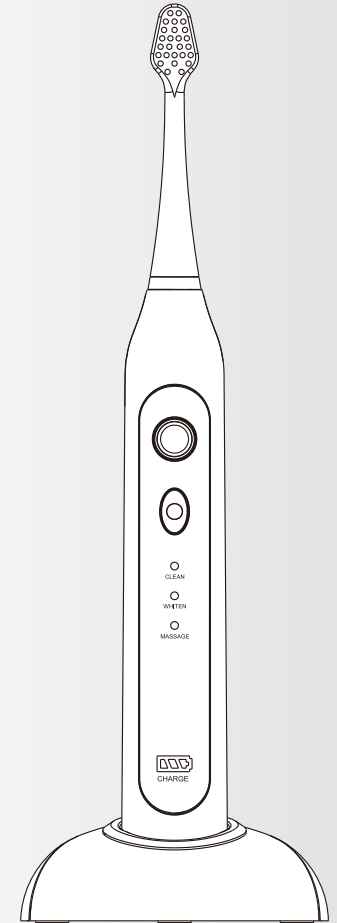ENG  
RO  
BG  
HU

**Dr. Mayer**  
**Life and Health**

E-mail: [office@doctor-mayer.com](mailto:office@doctor-mayer.com)

Web: [www.doctor-mayer.com](http://www.doctor-mayer.com)

Before operating this product, please read these instructions completely and save this manual for future use.

## CONTENTS

IMPORTANT SAFEGUARDS.....	1
DANGER-WATER.....	1
DANGER-VOLTAGE.....	1
DANGER-WARNINGS.....	1
MEDICAL-WARNINGS.....	3
Configuration.....	4
Getting Started.....	5
How To Use.....	6
Maintenance.....	7
Troubleshooting.....	8
Disposal.....	9

# IMPORTANT SAFEGUARDS

The Power toothbrush is safe and effective when used according to directions. Like any electrical appliance, however, basic safety precautions should be followed, especially when children are present.

READ ALL INSTRUCTIONS BEFORE USING:

## DANGER-WATER

To reduce the risk of burns, electrocution, fire, or injury to persons:

- ❶ Do not place or store the charger where it can fall or be pulled into a bathtub, sink or toilet.
- ❷ Do not place or drop the charger into water or other liquid.
- ❸ Do not reach for a product that has fallen into water or other liquid. Unplug immediately.
- ❹ Do not use while bathing.
- ❺ Unplug charger before cleaning and assure that it is dry before plugging it into an electrical outlet.

## DANGER-VOLTAGE

To reduce the risk of burns, electrocution, fire, or injury to persons:

- ❶ Do not plug the charger into an outlet with a voltage other than specified on the bottom of the charger.
- ❷ This product is designed to operate within a range of 100-240 volts.
- ❸ Voltage converters do not guarantee compatibility.
- ❹ Do not force the plug into an outlet.
- ❺ If the plug does not fit easily into the outlet, discontinue use.
- ❻ Never wrap the cord tightly around the appliance. Never allow the power supply cord to be pulled, twisted, or severely bent. Untwist the cord before use.

## DANGER-WARNINGS

To reduce the risk of burns, electrocution, fire, or injury to persons:

- ❶ This product is not intended for use by persons (including children) with

reduced physical, sensory or mental capabilities or lack of experience and knowledge, unless they have been given adequate supervision or instruction concerning use of the product by an adult responsible for their safety.

- ② Children should be supervised to ensure they do not play with the product.
- ③ Use this product only for its intended use as described in this manual. This product is designed to clean your teeth, gums and tongue only. Do not use any attachments or brush heads other than those recommended by the manufacturer.
- ④ The toothbrush is a personal care appliance and is not intended for use on multiple patients in a dental office or institution.
- ⑤ If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpastes), thoroughly clean your brush head and handle to avoid cracking the plastic.
- ⑥ Never operate this product if it has a damaged cord, plug or charger. Do not use this product if it is not working properly, if it has been dropped or damaged or dropped into water. Immediately discontinue use if the product appears damaged in any way.
- ⑦ Keep the product, including power cord away from heat surfaces.
- ⑧ If you experience discomfort or pain in your mouth, discontinue use of this product and contact a dentist.
- ⑨ Do not attempt to remove the battery unless the product life is end.
- ⑩ Do not clean the toothbrush in the dishwasher.
- ⑪ This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician.
- ⑫ Do not attempt to defeat this safety feature.
- ⑬ Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- ⑭ Do not leave appliance when plugged in. Unplug from outlet when not in use and before servicing.
- ⑮ Do not pull or carry by cord, use cord as a handle, close a door on cord, or pull cord around sharp edges or corners. Do not run appliance over cord.
- ⑯ Do not handle plug or appliance with wet hands.
- ⑰ Turn off all controls before unplugging.


## MEDICAL WARNINGS

- ① If you have had oral surgery or gum surgery within the past two months, consult your dental provider before using the product.
- ② Consult your dental provider if excessive bleeding occurs while using this product or minor bleeding continues to occur after one week of use.
- ③ Consult your physician or implant device manufacturer prior to using this product if you have questions concerning the combined use of this product with a pacemaker or other implanted medical device(s).
- ④ Consult your physician prior to using the toothbrush if you have any medical concerns.
- ⑤ Always unplug this product immediately after using.
- ⑥ Never use while sleeping or drowsy.


**SAVE THESE INSTRUCTIONS**

## Configuration


### Handle
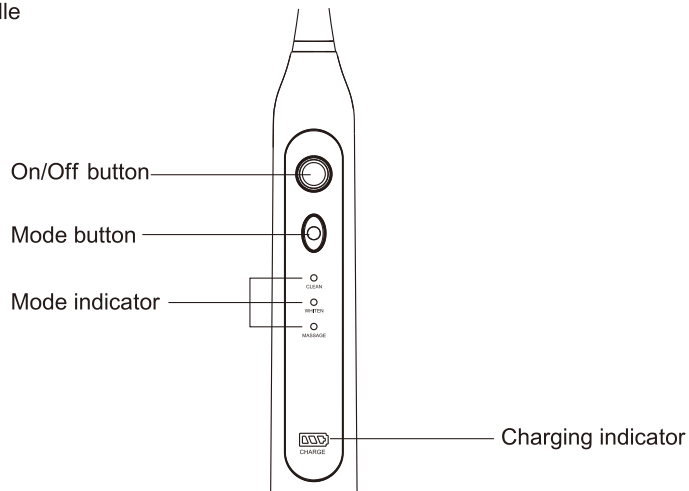
### Charger
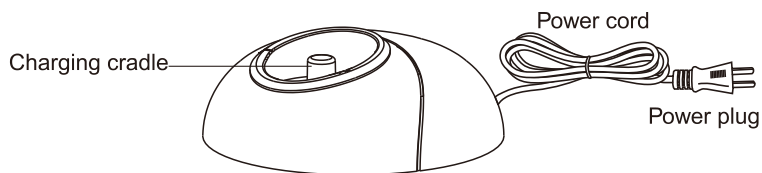
### Brush heads
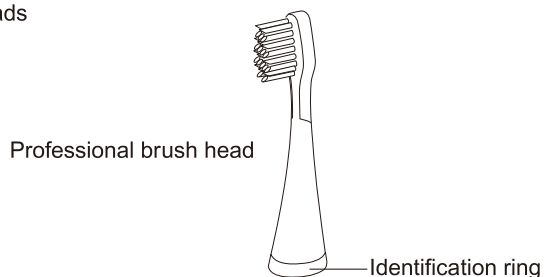
## Getting Started



### Charging the toothbrush

Be sure to turn off the toothbrush before charging.

Plug the power plug into the outlet. Place the handle in the charger. The charger indicator light will be lit, indicating the handle is charging. After the battery is fully charged, all the 3 sections of charge indicator light will turn solid green, indicating the toothbrush is fully charged.

### Notes:

- 1 It takes at least 24 hours to charge the battery fully. When fully charged, the toothbrush can be used 20 pre-timed brush cycles, or approximately 10 days, if brushing twice daily.
- 2 The toothbrush must be charged for 24 hours, if:
  - a. this is the first time using the toothbrush;
  - b. the toothbrush has not been used for more than three months; or
  - c. the charge is drained from use.
- 3 To keep the battery fully charged, we recommend that you keep the toothbrush in the charger when not in use.
- 4 During charging, the toothbrush will get warm, but this is normal.
- 5 The operating duration will become shorter as the battery

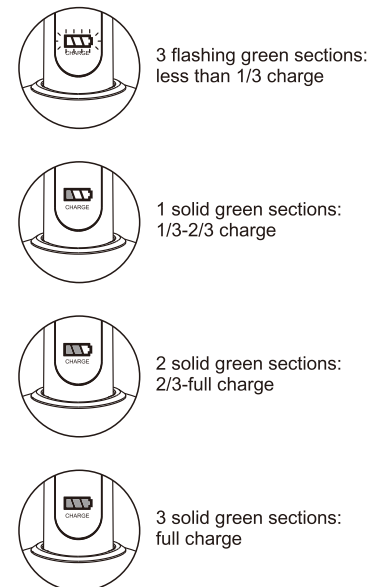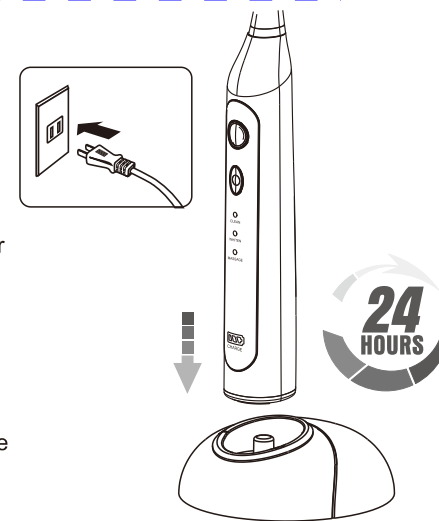gets older. If the toothbrush operates only for 1 short time after being fully charged, this means that the battery is nearing the end of its life cycle. Contact service center for support.

## How To Use

- ❶ Attach a brush head to the handle. There will be about 0.5mm space between the brush head and the handle. (Fig 1)
- ❷ Press the mode button to select mode. If you know the mode used last time, you need not select mode again. There are 3 modes.
  - a. CLEAN: Standard mode for superior teeth cleaning.
  - b. WHITEN: Teeth cleaning with polishing and whitening.
  - c. MASSAGE: Pulsates for gentle gum stimulation.
- ❸ Wet the bristles and apply a small amount of toothpaste.
- ❹ Place the bristles against the teeth at a slight angle towards the gum line. (Fig 2)
- ❺ Press the on/off button to turn on the toothbrush.
- ❻ Apply light pressure and gently move the brush head slowly across the teeth in a small back and forth motion.
- ❼ In order to maximize cleaning

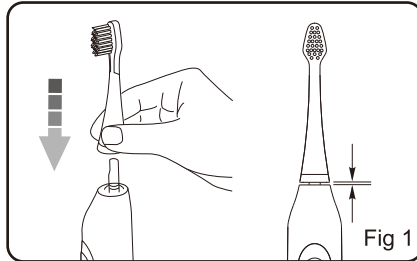
Fig 1

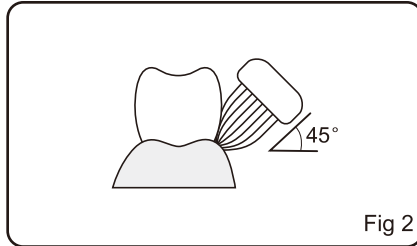
Fig 2

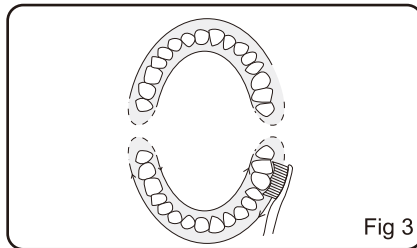
Fig 3

effectiveness, divide your mouth in 4 sections. Start with section 1 and brush each section for 30 seconds. The brushing action will pause between each 30 second interval, indicating when to move to the next section. The toothbrush will automatically stop after 2 minutes of use. (Fig 3)

- ❽ Take down the brush head from handle. (Fig 4)
- ❾ Rinse both the brush head and handle under running water and wipe them dry. (Fig 5)
- ❿ Put the handle back into the charger.

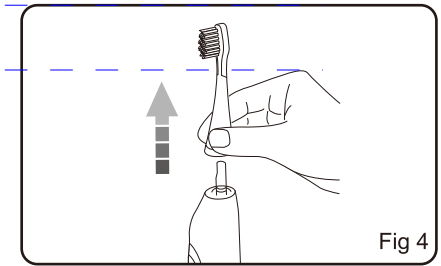
Fig 4

Fig 5

## Maintenance

Be sure to remove the plug from the electrical outlet before cleaning the product.

Charger

- ❶ Using a soft cloth, wipe off any dirt on the charger exterior surfaces.
- ❷ Clean the plug periodically especially the metal plate ensure it remains in good working order. The accumulation of dust or dirt on the plug can cause an electric shock, overheat or fire.

## Toothbrush

- ❶ Remove the brush head.
- ❷ Wipe all surface areas with a damp cloth and place the handle back into the charger. Do not submerge the handle into water.
- ❸ Rinse the brush head and bristles after each use.
- ❹ Do not clean in the dishwasher.

## Troubleshooting

Problem	Possible cause	Solution
The toothbrush doesn't work.	You have just purchased the product or it has not been used for more than 3 months.	Charge the handle for at least 24 hours.
Weak or diminished vibration.	The toothbrush is pressed hard against your teeth.	Place the brush head gently against your teeth.
The toothbrush operates only for a short period of time even immediately after being charged.	The battery has reached the end of its life cycle.	Contact service center for support.
	Charging time was insufficient.	Charge the handle for at least 24 hours.




## Disposal

The power source of this product is nickel-metal hydride battery which can be recycled. At the end of the product's life, always remove the battery before you discard the product, and recycle or dispose of the battery according to local waste management requirements. Do not throw away the appliance with normal household waste at the end its life. Instead, hand it in at an official waste collection point for recycling. By doing above all, you help to preserve and protect the environment. If you have recycling questions, contact your local waste management office.

### Removing and recycling the battery

Do not disassemble the product at any time other than when disposing of it. Doing so earlier will destroy the product. Observe basic safety precautions.

- ❶ Before removing the battery, be sure to completely drain the unit of power. Please turn on the device by press the on/off button, and let it work until it does not work.
- ❷ Use a screwdriver to remove the screw and battery door.
- ❸ Hold the handle upside down, press the toothbrush shaft, and remove the two hooks at the same time to release the internal components of the handle.
- ❹ Insert the screwdriver next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier.
- ❺ Recycle or dispose of the battery properly. Ask your local waste management office if they accept the remaining toothbrush components for recycling.

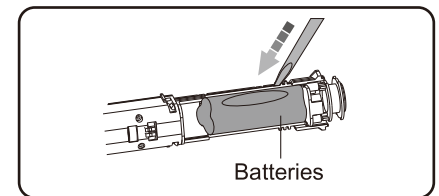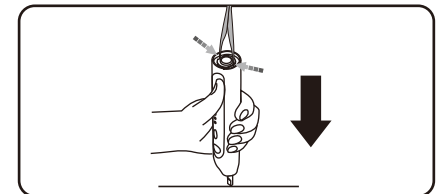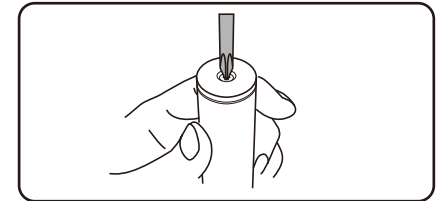Be careful not to short-circuit the negative and positive poles of the battery. Do not connect the appliance to the mains again after you have removed the battery.

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