
Countertop Water Flosser

Instruction Manual

WT 6000
Dr. Mayer
LIFE & HEALTH

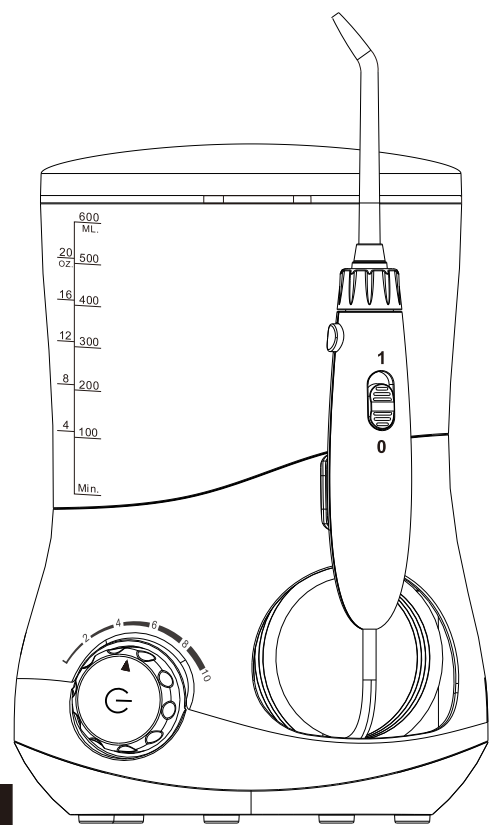ENG
RO
BG
HU

Include: Water flosser x1PC, Adaptor x1PC, Standard Nozzle x3PCS, Tongue Scraper x1PC, Instruction Manual x1PC

IMPORTANT SAFEGUARDS

Please read the safety instructions before using this product. The safety instructions provided here help you to use this product correctly, and also to prevent you from injuring yourself or others.

When using electrical products, especially when children are present, the basic safety precautions should always be followed.

READ ALL INSTRUCTIONS BEFORE USING

DANGERS

- Toddlers, infants, Children should not be allowed to use this device. Please keep it out of reach of children.
- Do not use if the power cord or plug was damaged.
- Always unplug the device after using.
- Do not touch plug with wet hand.
- Do not place the device in water or other liquids.
- Do not place the device where it may fall or be pulled into a tub or sink.
- Do not use the device while bathing.
- Do not touch the product which has fallen into water, Unplug immediately.


WARNINGS

- The device can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under supervision.
- Children should not play with the appliance.
- The appliance is only to be used with the power supply unit provided with the appliance.
- Keep the device away from heated surface.
- Do not use other tips expect recommend parts.
- Please do not use the device if you lost the nozzle.
- Keep nozzle away from small children. The nozzle is small enough to be swallowed and can be become lodged in the throat.
- Do not direct water under the tongue, into the ear, nose or other delicate areas.
- Do not use for any other purpose except oral cleaning.
- Do not fill the tank with warm water over than 40 °C.
- Certain types of mouthwash may damage the device by causing the water tank to crack. Therefore, do not use mouthwash in the water tank.
- People who experience tooth or gum pain may not be able to use this device and should consult their dentist.
- Do not press the nozzle too close against your teeth or gums, doing so may injure your gums.
- For household use only, not for outside use.
- Please fill the tank with pure water.

SAVE THESE INSTRUCTIONS


PARTS IDENTIFICATION

Water Flosser

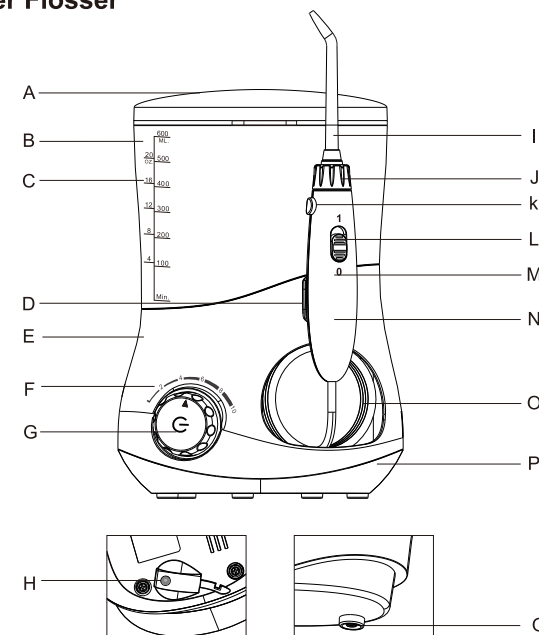- | | |
|---|-------------------|
| A: Nozzle Case | I: Nozzle |
| B: Water Tank | J: Nozzle knob |
| C: Volume Mark | K: Release Button |
| D: Handle Holder | L: Pause Switch |
| E: Unit Base | M: On/off Mark |
| F: Pressure Mark | N: Handle |
| G: Main Button | O: Pipe |
| (Pressing: On/Off control Rotation: Pressure adjust) | P: Unit Bottom |
| H: DC Jack Plug | Q: Tank Valve |

Accessories

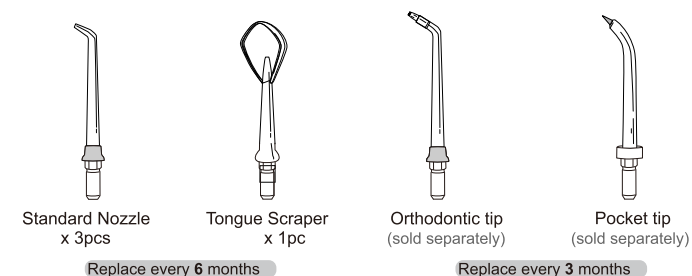HOW TO USE

GETTING STARTED

Place the device on a flat surface without any sharp objects and make sure the cord is dry and has no damage. The device is powered by a special adaptor with DC jack connector. (Input: 100-240V~, 50/60Hz; Output: 12V , 1.5A).

1. Insert the plug into the DC jack on the bottom of base and plug the adaptor cord into an electrical outlet. (Fig.1)
2. Take the Nozzle Case off.
3. Take the Water Tank off to fill with proper liquid, no warmer than 40 °C, do not pour in salt, tooth-paste, or chemicals, then replace the tank onto the top of main unit and press down firmly. (Fig.2-1&Fig2-2)

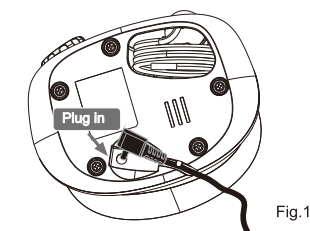
Fig.1

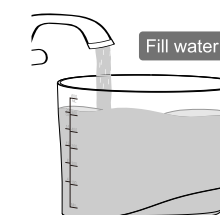
Fig.2-1

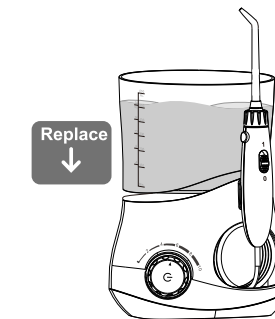
Fig.2-2

Inserting & Removing Nozzle

Inserting Nozzle

Insert tip into the center of the handle, audibly indicating when the nozzle attached to the handle in place. And also try to pull the tip out gently to confirm the tip is locked firmly. (Fig.3)

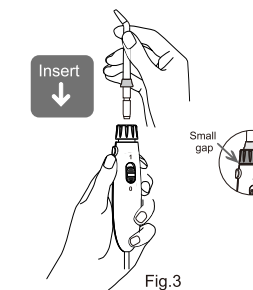
Fig.3

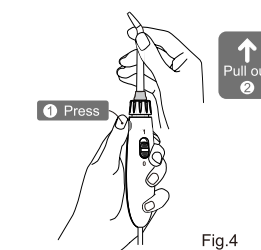
Fig.4

Removing Nozzle

Press the Release Button on the top end of the handle, then pull the nozzle out lightly by hand. (Fig.4)

How to use the water flosser

Note:

Do not power on the device before placing the nozzle in your mouth. First use, we suggest you choose the lower pressure, gradually increase the water pressure over time to the setting you need.

1. Hold the handle and point the nozzle to sink, then press down the Main Button to power the device on. (Fig.6)
2. Rotate the Main Button to adjust water pressure, after setting, push the Pause Switch to "0" position to pause the water flow. (Fig.7-1&Fig.7-2)

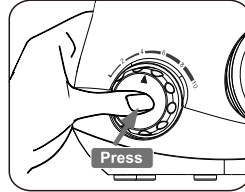
Fig.6

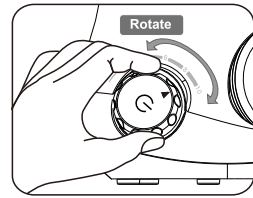
Fig.7-1

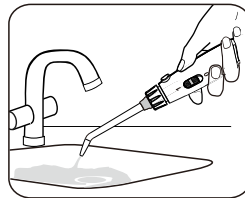
Fig.7-2

3. Adjust the direction of nozzle by rotating the nozzle knob;
4. Bend over to the sink and place the tip in your mouth, and then push the Pause Switch to "1" position to enjoy irrigating. (Slightly close your mouth to avoid the water splashing.) (Fig.8)

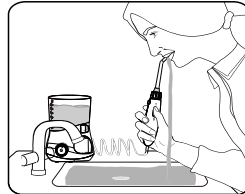
Fig.8

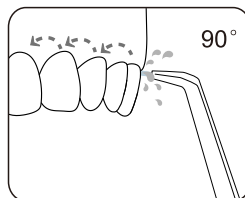
Fig.9

5. Position the water jet flow at 90-degree angle to the teeth and gums. (Fig.9)

6. Move the water jet slowly along the teeth. The water flow will short pause at **30 seconds**, which reminds you to clean the other arch of your mouth and another short pause at **1-minute** to remind you the recommended usage time is up. The unit will continue to work till **2-minute** cycle is complete or you turn it off. (Fig.10)

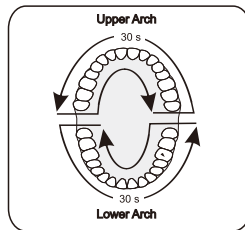
Fig.10

How to end use

You can end irrigating through the below methods:

- 1: Auto shut-off turns off after **2-minute** cycle is complete;
- 2: Turn it off by pressing the Main Button;
- 3: Unplug the adaptor directly when in any emergency;

After use

1. Take off the water tank and dispose of water in it.
2. Push the Pause Switch to "1" position, and power on to drain the water inside the device, then push the Pause Switch to "0" position again.
3. Turn off the device and wipe it with a soft dry cloth. (Fig.11)
4. Press down the Release Button, and remove the nozzle, then place it into the nozzle case. (Fig.12)
5. Unplug the adaptor.

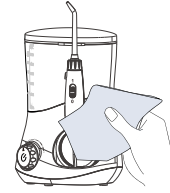
Fig.11

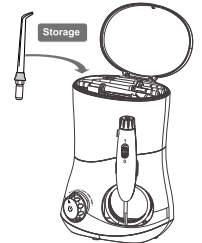
Fig.12

FUNCTIONAL TIP

Tongue Scraper

When you want to clean your tongue, please gently place the tongue scraper onto the center midline of your tongue surface about halfway back with light water pressure setting first, and gradually increase the pressure as needed. (Fig.13)

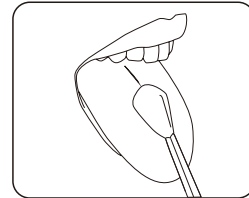
Fig.13

Orthodontic Tip (sold separately)

The orthodontic tip is specially designed for people wearing braces: it will clean the hard to reach areas around and in braces.

- 1: Aim water stream at 90 degree to the tooth at the gum-line, start on the last molar (back tooth);
- 2: Slowly slide the tip along with the gum-line, and pausing briefly to use the brush bristles to clean the area between teeth and all around orthodontic braces.
- 3: Repeat until you have finished cleaning the inside and outside of all your teeth. (Fig.14)

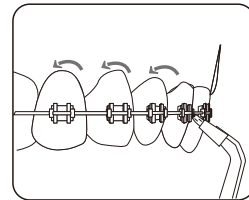
Fig.14

Pocket Tip (sold separately)

The pocket tip is specifically designed for the cleaning of periodontal pockets. To use the pocket tip, it is a MUST to set the device to the lowest pressure setting.

1. Place the soft tip against a tooth at 45-degree angle and gently place the tip under the gum-line, then insert into the pocket. (Fig.15)

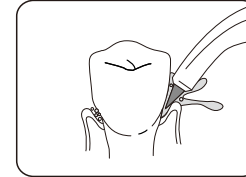
Fig.15

2. Turn the unit on and continue tracing along the gum-line until all periodontal pockets are treated. (Fig.16)

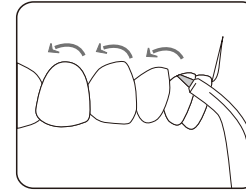
Fig.16

Note:

This tip must be used only with the lowest water pressure to extend life of soft tip as it is easily broken.

MAINTENANCE

Use only water or neutral cleaning agent to clean the device. Never use caustic or abrasive agents (e.g., vinegar or lime scale removers), since they may damage the device. Do not wash using hot water of over 40 °C. Keep the device away from high temperature area and avoid direct sunlight.

Main unit

1. Remove the water tank and drain the remains water inside the device.
2. Clean by using a soft cloth with liquid soap and dry with a soft dry cloth.
3. Do not submerge the main unit in water to clean.

Handle

1. Wash with water and dry up with a soft cloth.
2. Do not bend, pull or twist the water hose.

Water tank

1. Wash with water and wipe with a soft dry cloth.
2. Make sure the water tank is inverted and firmly fixed to store.

Nozzle

1. Wash with water and wipe with a soft dry cloth.
2. Place the nozzle in the storage compartment.
3. Replace once every 6 months for Standard nozzle and tongue scraper, and once every 3 months for other tips.

Charging adaptor

Wipe the metal part of the power plug with a dry cloth once every 6 months.

TROUBLESHOOTING

| Problem | Possible Cause | Solution |
|----------------------------|---|--|
| Unit does not work | Motor not running | Make sure the adapter is firmly connected to the outlet and the unit. Make sure the outlet and unit power is on |
| | Dirt stuck in the nozzle | Clean or replace with a new nozzle |
| Water pressure is very low | The nozzle is deformed | Replace with a new nozzle |
| | Dirt stuck in the tank valve | Clean the valve and fill the tank with pure water |
| Water does not come out | The tank is empty | Fill the tank up |
| | Water pressure under the lowest setting | Slide the switch to adjust pressure |
| | Tank valve not open | Press down the tank, or reassemble the tank valve |
| Tank leaks | The valve was lost | Replace with a new tank valve |
| | The valve isn't install in place | Take off the valve and reinstall again |
| Tear in Pocket tip head | Tip head is worn | Replace with a new nozzle |

DISPOSAL

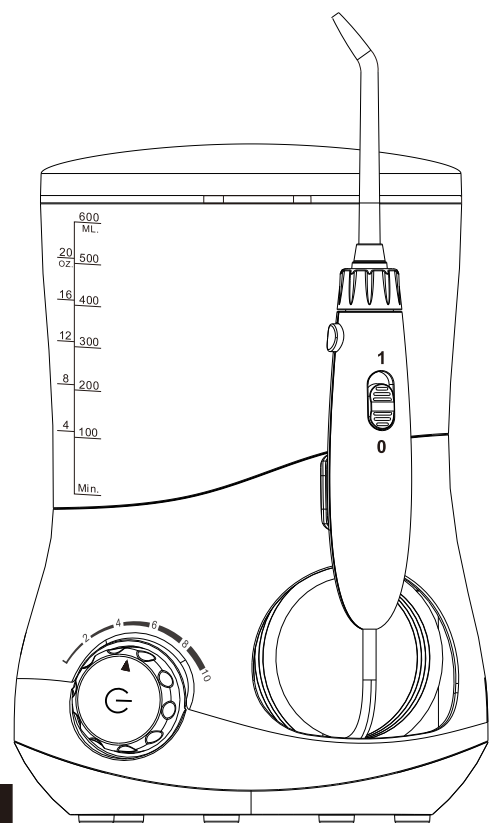
Do not disassemble the device at any time other than when disposing of it. Doing so earlier will destroy the product. The appliance must be disconnected from the power supply when disassembling the device.

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